

## A Protocol to Measure the Cytoplasmic Calcium in *Arabidopsis* Guard Cells

Li Li, Feng Lin, Yana Qu and Qun Zhang\*

College of Life Sciences, State Key Laboratory of Crop Genetics and Germplasm Enhancement, Nanjing Agricultural University, Nanjing, China

\*For correspondence: [zhangqun@njau.edu.cn](mailto:zhangqun@njau.edu.cn)

**[Abstract]** Cytoplasmic calcium ( $[Ca^{2+}]_{cyt}$ ) acts as a stimulus-induced second messenger in multiple signal transduction cascades (Allen *et al.*, 1999). In plant cells, a dramatic and readily assayed response to stimulus is the change of stomatal aperture. Changes in  $[Ca^{2+}]_{cyt}$  of stomatal guard cells were involved in stomatal movement in response to various stimuli and cellular processes. In general, there are two available ways to measure  $[Ca^{2+}]_{cyt}$  in guard cells, i.e., loading of calcium-sensitive fluorescence dyes such as fluo-3 AM and fura-2 or expressing genetically encoded calcium indicators such as yellow cameleon (Krebs *et al.*, 2012). In this protocol, we aim at describing the experimental procedure to record  $[Ca^{2+}]_{cyt}$  fluctuation in guard cells with loading of fluo-3 AM upon ABA or PA treatment combining with fluorescence imaging performed with confocal laser scanning microscope.

### Materials and Reagents

1. Plant materials: *Arabidopsis thaliana* ecotype Columbia-0 (wild type) was obtained from ABRC at Ohio State University (Columbus)
2. Glycine,N-[4-[6-[(acetyloxy)methoxy]-2,7-dichloro-3-oxo-3H-xanthen-9-yl]-2-[2-[2-[bis[2-[(acetyloxy)methoxy]-2-oxyethyl]amino]-5-methylphenoxy]ethoxy]phenyl]-N-[2-[(acetyloxy)methoxy]-2-oxyethyl]-, (acetyloxy) methyl ester (Fluo-3, AM) (Life Technologies, Invitrogen™, catalog number: F-14218)
3. Potassium chloride (KCl) (Sangon Biotech, catalog number: PB0440)
4. Calcium chloride dihydrate ( $CaCl_2 \cdot 2H_2O$ ) (Sangon Biotech, catalog number: C0556)
5. 2-(N-morpholino) ethanesulfonic acid (MES) (free acid, monohydrate) (Sigma-Aldrich, catalog number: M3671)
6. Potassium hydroxide (KOH) (Sangon Biotech, catalog number: PT1159)
7. Egtazic acid, Glycol ether diamine tetraacetic acid (EGTA) (Sangon Biotech, catalog number: ED007)
8. Chloroform (Sinopharm Chemical Reagent, catalog number: 10006818)
9. Sodium hypochlorite solution (Sangon Biotech, catalog number: S1944)
10. Agar (Sangon Biotech, catalog number: AJ637)
11. Murashige & Skoog (MS) Basal Medium w/Vitamins (*PhytoTechnology Laboratories*®, catalog number: M519)
12. Abscisic acid (ABA) (Sigma-Aldrich, catalog number: A1049)

13. Phosphatidic acid (PA) (Avanti Polar Lipid, catalog number: 840858C)
14. Dimethyl sulfoxide (DMSO) (Sigma-Aldrich, catalog number: D2650)
15. Ethanol anhydrous (Sangon Biotech, catalog number: ET0737)
16. Tris (hydroxymethyl) aminomethane (Tris) (Sigma-Aldrich, catalog number: 252859)
17. Epidermal buffer (see Recipes)
18. Fluo-3 AM stock solution (see Recipes)
19. Loading buffer (see Recipes)
20. ABA stock solution (see Recipes)
21. ABA working solution (see Recipes)
22. PA working solution (see Recipes)

### **Equipment**

1. Graduated centrifuge tubes (1.5 ml)
2. Petri dish (6 cm diameter)
3. Cover glass (20 mm \* 20 mm, 0.13-0.17  $\mu\text{m}$  thickness)
4. Slides (76 mm \* 26 mm, 1-2 mm thickness)
5. Forceps
6. Pipettor and matched tips
7. Growth chamber (Percival, model: I-41LL)
8. Zeiss LSM780 confocal laser scanning microscope
9. Microprocessor pH meter (HANNA Instruments, model: HI-2211)
10. Ultrasonic cell crusher (Scientz, model: JY92 II)

### **Software**

1. Zeiss confocal software ZEN 2012 (blue edition)
2. Microsoft Excel

### **Procedure**

#### **A. Sample preparation**

1. Seeds of *Arabidopsis thaliana* were surface-sterilized by soaking in 75% ethanol for 3 min and 2% sodium hypochlorite for 10 min. The seeds were washed for 3 to 5 times with sterile distilled water, and then were placed on MS medium (1% agar). *Arabidopsis* seedlings were grown in a growth chamber at 160  $\mu\text{mol}/\text{m}^2/\text{s}$  light intensity and 14/10 h (24/21 °C) day/night regimes. The cotyledons of 7-day-old seedlings were cut from the plants and used for experiments.
2. The cotyledons were incubated with epidermal buffer at 24 °C for 2 h in growth chamber (160  $\mu\text{mol}/\text{m}^2/\text{s}$  light intensity) to induce stomatal opening prior to various

treatments.


3. The samples were then incubated in 10  $\mu\text{M}$  fluo-3 AM loading buffer (diluted from 5 mM stock solution in 10 mM MES-KOH pH 6.1) for 2 h at 4  $^{\circ}\text{C}$  in darkness. Higher concentration of fluo-3 AM at room temperature for incubation may result in false positive signals from other ions, for example  $\text{Mn}^{2+}$ ,  $\text{Zn}^{2+}$  and  $\text{Pb}^{2+}$ .
4. The samples were rinsed with loading buffer for 3 times to remove fluo-3 AM and were kept in growth chamber (24  $^{\circ}\text{C}$ ) for 1 h before ABA or PA addition.
5. For ABA or PA treatment, the cotyledons were incubated with 10  $\mu\text{M}$  ABA or 50  $\mu\text{M}$  PA working solution, and observed by the confocal laser scanning microscope.

#### B. $[\text{Ca}^{2+}]_{\text{cyt}}$ imaging

1.  $[\text{Ca}^{2+}]_{\text{cyt}}$  imaging was performed with a 488 nm excitation and a 510-545nm emission filter to record  $[\text{Ca}^{2+}]_{\text{cyt}}$  images. Increases in  $[\text{Ca}^{2+}]_{\text{cyt}}$  occurred after ABA or PA addition. Fluorescence intensity of Fluo-3 dye increases when bound to  $[\text{Ca}^{2+}]$  in the cytosol. The  $K_d$  of Fluo-3 for calcium is 325 nM. Fluo-3 can also detect calcium near the chloroplasts membrane and plasma membrane, and display higher fluorescence in cytosol because of higher concentration of calcium. Changes in  $[\text{Ca}^{2+}]_{\text{cyt}}$  were visualized in pseudocolour as indicated by the color bars. The color bar showed the pixel intensity as the arrow indicated (Figure 1).
2. The imaging parameters were as follows: Plan-Apochromat 10 $\times$ /0.45 M27 objective, image dimension 1,024 x 1,024, pinhole 3 Airy unit, scanning speed 1936.28 ms. The images were acquired at the time indicated. Fluorescence pixel intensity values in guard cells were measured using the Zeiss confocal software ZEN 2012 (blue edition). The pixel intensity value of background is about 4.5 for subtraction.

*Please note that the laser power, gain, offset, zoom, scan speed, etc. should not be changed across experimental conditions if using intensity to distinguish phenotype. Set values at the brightest sample and continue imaging with same conditions to get an appropriate dataset. If the fluorescence is high enough, it can be performed with a fluorescence scope. The confocal scanning will be more sensitive and reproducible.*

**Representative data**

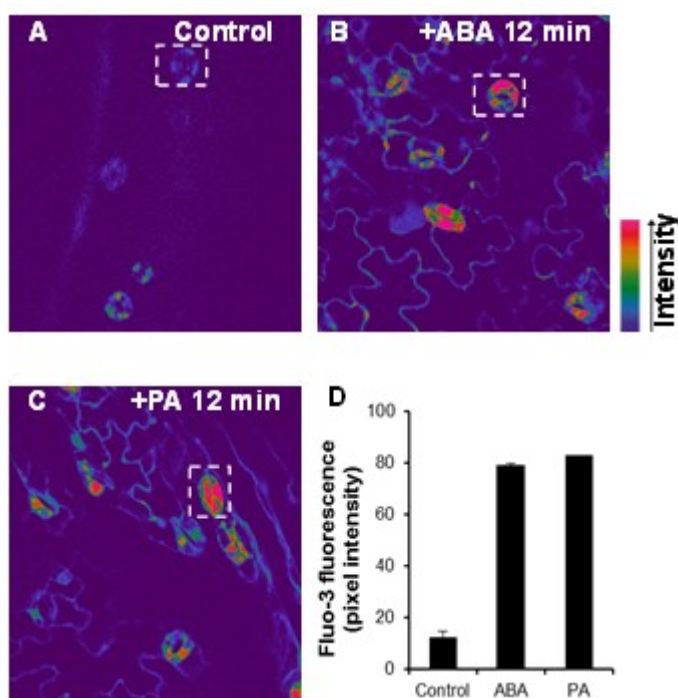

**Figure 1. Elevation of cytoplasmic calcium induced by ABA and PA in *Arabidopsis* guard cells.** The confocal images of  $[Ca^{2+}]_{cyt}$  in guard cells were monitored by fluo-3 AM dye before (A) or after treatment with 10  $\mu$ M ABA (B) or 50  $\mu$ M PA (C) for 12 min. The color bars showed the intensity as the arrow indicated. D. Changes in the relative levels of  $[Ca^{2+}]_{cyt}$  in ABA- or PA-treated guard cells. Regions of interests used to measure the intensities were indicated by white rectangle with background subtraction. Values are the mean  $\pm$  SD ( $n = 50-60$  from not less than 10 cotyledons) from three independent experiments.

**Recipes**

1. Epidermal buffer
  - 10 mM KCl
  - 0.2 mM  $CaCl_2$
  - 0.1 mM EGTA
  - 10 mM MES-KOH
  - Adjust the pH to 6.15 with KOH (0.1 M)
2. Fluo-3 AM stock solution
  - 5 mM Fluo-3 AM in DMSO
3. Fluo-3 AM Loading buffer
  - 10 mM MES
  - Adjust the pH to 6.1 with KOH (0.1 M)
  - 10  $\mu$ M Fluo-3 AM

4. ABA stock solution  
Dissolve 26.432 mg ABA powder in 1 ml ethanol to a concentration of 100 mM
5. ABA working solution (for treatment)  
Add ABA stock solution (10,000x) to the epidermal buffer to a concentration of 10  $\mu$ M ABA
6. PA working solution (for treatment)  
PA in chloroform were dried under N<sub>2</sub> and suspended in 1 ml of epidermal solution by sonication for preparing the stock solution (50 mM, 1,000x). The PA stock solution was diluted into 50  $\mu$ M with epidermal solution before treatment.

### **Acknowledgments**

The methods were adapted from (Jiang *et al.*, 2014). Techniques were also adapted from all of the references cited. This work was supported by grants from National Basic Research Program of China (31100194 and 31470364) and the Fundamental Research Funds for the Central Universities (KYZ201423) to Q Zhang.

### **References**

1. Allen, G. J., Kwak, J. M., Chu, S. P., Llopis, J., Tsien, R. Y., Harper, J. F. and Schroeder, J. I. (1999). [Cameleon calcium indicator reports cytoplasmic calcium dynamics in \*Arabidopsis\* guard cells.](#) *Plant J* 19(6): 735-747.
2. Jiang, Y., Wu, K., Lin, F., Qu, Y., Liu, X. and Zhang, Q. (2014). [Phosphatidic acid integrates calcium signaling and microtubule dynamics into regulating ABA-induced stomatal closure in \*Arabidopsis\*.](#) *Planta* 239(3): 565-575.
3. Krebs, M., Held, K., Binder, A., Hashimoto, K., Den Herder, G., Parniske, M., Kudla, J. and Schumacher, K. (2012). [FRET-based genetically encoded sensors allow high-resolution live cell imaging of Ca\(2\)\(+\) dynamics.](#) *Plant J* 69(1): 181-192.