

## Spot Assays for Viability Analysis of *Cyanobacteria*

Anja Katharina Dörrich<sup>1</sup> and Annegret Wilde<sup>2,3\*</sup>

<sup>1</sup>Institute for Microbiology and Molecular Biology, Justus-Liebig-University, Giessen, Germany;

<sup>2</sup>Institute of Biology III, University of Freiburg, Freiburg, Germany; <sup>3</sup>Centre for Biological Signalling Studies (BIOSS), University of Freiburg, Freiburg, Germany

\*For correspondence: [annegret.wilde@biologie.uni-freiburg.de](mailto:annegret.wilde@biologie.uni-freiburg.de)

**[Abstract]** *Cyanobacteria* are prokaryotic organisms performing oxygenic photosynthesis. The cyanobacterium *Synechocystis* sp. PCC 6803 is a model organism for the study of photosynthesis, gene regulation and biotechnological applications because it is easy to manipulate genetically. Moreover, this cyanobacterium can grow photoautotrophically as well as chemoheterotrophically in the dark utilizing glucose. Microbiologists often use optical density measured with a spectrophotometer for the comparison of growth performance of different strains in liquid cultures. Because *Synechocystis* sp. PCC 6803 (especially motile strains) tend to form aggregates under stress conditions this method might be not suitable for evaluation of different strains under different growth conditions. In addition, many labs are not well equipped with standardized photobioreactors and illumination facilities to ensure reproducibility of growth curves. Here, we describe a highly reproducible spot assay for viability analysis of cyanobacterial strains.

### **Materials and Reagents**

1. *Synechocystis* sp. PCC 6803 strain (wild type obtained from S. Shestakov, Moscow State University, Russia)
2. CaCl<sub>2</sub>·2H<sub>2</sub>O (Carl Roth, catalog number: 5239.1)
3. Citric acid (Carl Roth, catalog number: X863.2)
4. NaNO<sub>3</sub> (Carl Roth, catalog number: A136.2)
5. MgSO<sub>4</sub>·7H<sub>2</sub>O (Carl Roth, catalog number: P027.1)
6. Na<sub>2</sub>-EDTA (Carl Roth, catalog number: 8043.2)
7. H<sub>3</sub>BO<sub>3</sub> (Carl Roth, catalog number: 6943.1)
8. MnCl<sub>2</sub>·4H<sub>2</sub>O (Carl Roth, catalog number: T881.3)
9. ZnSO<sub>4</sub>·7H<sub>2</sub>O (Carl Roth, catalog number: T884.1)
10. Na<sub>2</sub>MoO<sub>4</sub>·2H<sub>2</sub>O (Carl Roth, catalog number: 0274.1)
11. Co(NO<sub>3</sub>)<sub>2</sub>·6H<sub>2</sub>O (Carl Roth, catalog number: HN16.1)
12. CuSO<sub>4</sub>·5H<sub>2</sub>O (Carl Roth, catalog number: P024.1)
13. Na<sub>2</sub>CO<sub>3</sub> (Carl Roth, catalog number: P028.2)
14. Ammonium iron (III) citrate (Carl Roth, catalog number: P027.1)
15. K<sub>2</sub>HPO<sub>4</sub>·3H<sub>2</sub>O (Merck KGaA, Calbiochem<sup>®</sup>, catalog number: 105099)

16. 2-[(2-Hydroxy-1, 1-bis(hydroxymethyl)ethyl)amino]ethanesulfonic acid (TES) (Carl Roth, catalog number: 9137.3)
17. Bacto Agar (BD Bioscience, catalog number: 214010)
18. Syringes (50 ml volume)
19. Sterile syringe filters (Sarstedt Filtropur S 0.2  $\mu\text{m}$ , catalog number: 83.1826.001 and 0.45  $\mu\text{m}$ , catalog number: 83.1826)
20. Sterile 96-well microliter plates (*i.e.* VWR International, catalog number: 734-2328)
21. 100x BG11 medium (see Recipes)
22. Trace metal mix (see Recipes)
23. Stock Solutions for BG-11 medium (see Recipes)
24. 1x BG11 (2x BG11) (see Recipes)
25. BG11 agar (see Recipes)
26.  $\text{Na}_2\text{S}_2\text{O}_3$  solution (see Recipes)
27. 0.25 M  $\text{Na}_2$ -EDTA (see Recipes)

### **Equipment**

1. Wide neck Erlenmeyer flasks with 100 ml volume (Carl Roth, catalog number: C146.1)
2. Light source (18 W/840) (Philips, model: MASTER TL-D Super 80)
3. Laboratory shaker (one-dimensional orbital motion) (Heidolph Instruments GmbH, catalog number: 036130180)
4. UV/Vis spectrophotometer (Shimadzu, model: UV-2401PC)
5. Multichannel pipettes (8-channel pipette, 0.5-10  $\mu\text{l}$  and 10-100  $\mu\text{l}$ ) (Starlab ErgoOne<sup>®</sup>, catalog number: S7108-0510 and S7108-1100)
6. Petri dishes with square shape (120 x 120 x 17 mm) (Greiner Bio-One GmbH, catalog number: 688102)
7. Transmitted light scanner with adjustable cover (*i.e.* Epson perfection v700 photo, a good camera can be used as an alternative)
8. Laminar flow hood (*i.e.* Heraeus HeraSafe HS)
9. LI-190R Quantum Sensor

### **Procedure**

Every experimental step involving the handling of cyanobacterial cell cultures must be performed under a Laminar flow hood to prevent contamination!


1. Cyanobacterial liquid cultures were inoculated from agar plates in BG11 medium (Rippka *et al.*, 1979) and incubated with continuous shaking (150 rpm) at 30 °C under constant illumination with white light of 50  $\mu\text{mol photons m}^{-2} \text{s}^{-1}$  measured at the position of the agar plates for 2-3 days (measured with a LI-190R Quantum Sensor).

*Note: Even though some mutants carry antibiotic resistance gene cassettes inserted in their genome, cultures for the spot assays were always grown in medium lacking antibiotics to exclude effects from growth on selective media. There is no specific volume or OD needed, but the pre-cultures should not reach late stationary phase and their OD<sub>750</sub> should be similar.*

- The well-grown pre-cultures were diluted to OD<sub>750</sub> 0.2 in a final volume of 20 ml BG11 medium and were incubated further under the above mentioned conditions until they reached logarithmic phase (approximately OD<sub>750</sub> 0.6-0.7).

*Note: It is important to measure OD at a wavelength where no pigments absorb.*

- The logarithmically growing cultures were again diluted to OD<sub>750</sub> 0.2. With the freshly diluted cultures, a dilution series was performed with a multichannel pipette in a sterile microliter plate ranging from 10<sup>0</sup> to 10<sup>-7</sup> (Figure 1).

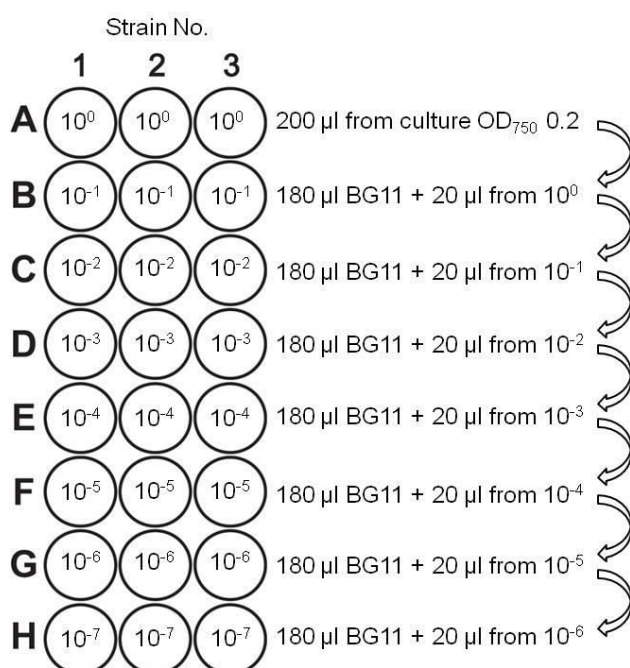


**Figure 1. Dilution series for viability analysis.** Rows B to H were pre-filled with 180 µl BG11 medium. After every dilution step, the cell suspension was thoroughly mixed with the pipette before transferring part of it into the next well.

- From each of the dilution steps, 5 µl was spotted in triplicates using the multichannel pipette onto square BG11 agar plates supplemented with 0.3% (w/v) Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub>.

*Note: The exact composition of the agar can be customized. If the experimenter wishes to observe viability under mixotrophic conditions, the BG11 agar can be i.e. supplemented with glucose. It is also possible to exclude iron or nitrogen from the media stock solutions. However do not forget to add Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub> to the plates. This is important for proper growth of the colonies (Thiel et al., 1989).*

5. The drops of liquid cultures were allowed to dry for approximately 15-20 min (closed plate, at room temperature on the table top). Afterwards, plates were incubated at 30 °C under constant illumination with white light at 50  $\mu\text{mol photons m}^{-2} \text{s}^{-1}$  for approximately 5 days or until no more colonies appeared in the highest dilution step.
 

*Note: To prevent phototactic movement of motile strains, the plates should be placed below the light source and lateral illumination from only one side should be avoided. In order to prevent dripping (in case the drops are not completely dry), the plates were incubated with the lid facing upwards for the first day. After the first day, plates were turned upside down to avoid evaporated water dripping on the spots. For investigation of light-dark related phenotypes, constant illumination can also be changed to light-dark-cycles. In our case, we investigated mutants of circadian clock associated genes and therefore chose a 12:12 h light:dark incubation. In this case, the starting  $OD_{750}$  for the dilution series was adjusted to 0.4 instead of 0.2 (step 2).*
6. Record the result (Figure 2) by scanning the plate or by taking a picture. In the best case, it may even be possible to count the single colonies within several spots (you may need a Binocular microscope).

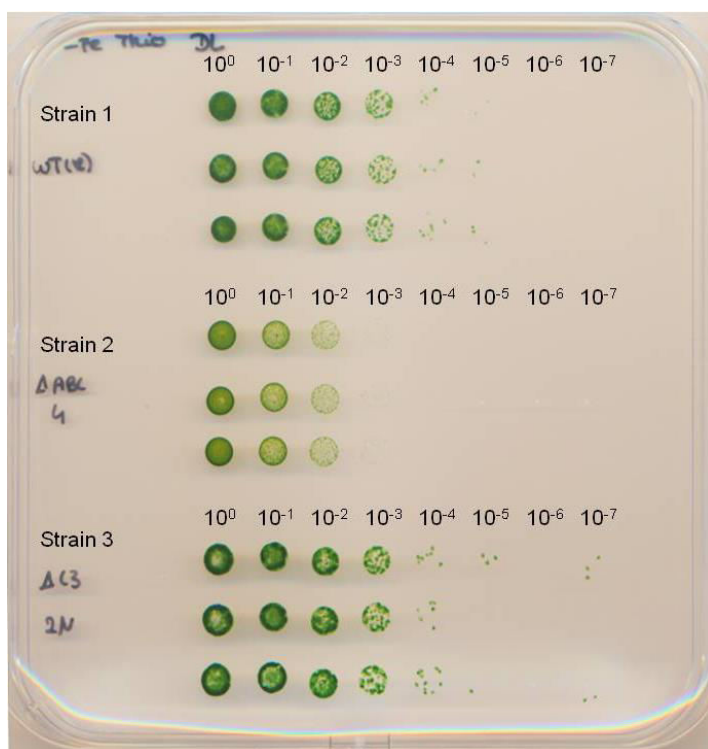

**Figure 2. Spot assay.** The dilution series of the cell cultures were spotted in triplicates onto square BG11 agar plates. The plate was incubated until no more colonies appeared in the highest dilution step. Strains have been published by Dörrich *et al.* (2014).

Further results generated by using this method were published by Dörrich *et al.* (2014).

## Notes

The spot assays described above are highly reproducible. In contrast to other methods like measuring the optical density of liquid cultures, which is often error-prone due to *i.e.* the formation of aggregates, we always obtained reliable and reproducible results. Within our group, several researchers have performed this method in two different laboratories and we did not detect any variations in the results. For reproducibility it is highly recommended to use Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub> in the plates (as described in Recipes) and to record the results after the same time. Composition of agar plates and light conditions can be customized to fit the experimenters' individual needs and may require optimization for different mutant strains.

## Recipes

### 1. 100x BG11 medium

CaCl <sub>2</sub> ·2H <sub>2</sub> O	3.6 g/L
Citric acid	0.6 g/L
NaNO <sub>3</sub>	149.58 g/L
MgSO <sub>4</sub> ·7H <sub>2</sub> O	7.49 g/L
0.25 M Na <sub>2</sub> -EDTA (pH 8)	0.56 ml/L

Prepare the solution in ddH<sub>2</sub>O and sterilize by autoclaving

### 2. Trace metal mix

H <sub>3</sub> BO <sub>3</sub>	2.86 g/L
MnCl <sub>2</sub> ·4H <sub>2</sub> O	1.81 g/L
ZnSO <sub>4</sub> ·7H <sub>2</sub> O	0.222 g/L
Na <sub>2</sub> MoO <sub>4</sub> ·2H <sub>2</sub> O	0.390 g/L
Co(NO <sub>3</sub> ) <sub>2</sub> ·6H <sub>2</sub> O	0.049 g/L

Prepare these solutions in ddH<sub>2</sub>O and sterilize by filtering using a syringe and sterile filter

Sterile aliquots can be stored at 4 °C

### 3. Stock Solutions for BG-11 medium

Na <sub>2</sub> CO <sub>3</sub>	20 mg/ml
Ammonium iron (III) citrate	6 mg/ml
K <sub>2</sub> HPO <sub>4</sub> ·3H <sub>2</sub> O	30 mg/ml
TES buffer (pH 8.0 adjust with NaOH)	1 M

Prepare these solutions in ddH<sub>2</sub>O and sterilize by filtering using a syringe and sterile filter

Sterile aliquots can be stored at 4 °C

### 4. 1x BG11 (2x BG11)

100 x BG11	10 ml/L (20 ml/L)
------------	-------------------

Na <sub>2</sub> CO <sub>3</sub> stock solution	1 ml/L (2 ml/L)
K <sub>2</sub> HPO <sub>4</sub> ·3H <sub>2</sub> O stock solution	1 ml/L (2 ml/L)
Ammonium iron (III) citrate stock solution	1 ml/L (2 ml/L)
Trace metal mix stock solution	1 ml/L (2 ml/L)
TES buffer (pH 8.0) stock solution	10 ml/L (20 ml/L)

Add sterile solutions to ddH<sub>2</sub>O and adjust the volume to 1 L

If you use non-sterile solutions and autoclave the medium afterwards there will be slight precipitation of salts. Before using the medium you should mix it thoroughly. This is recommended, when you want to avoid contaminations in your cultures, especially for long-term strain maintenance.

#### 5. BG11 agar

Prepare 1.5% bactoagar in ddH<sub>2</sub>O and autoclave.

The molten 1.5% bactoagar (approx. 80 °C) is mixed with the same volume of sterile 2x BG11 medium.

#### 6. Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub> solution

30% w/v Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub>

30 g Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub>

Add ddH<sub>2</sub>O to 100 ml

Filter sterilize (0.2 µm)

Stored at 4 °C

#### 7. 0.25 M Na<sub>2</sub>-EDTA (pH 8)

9.306 g Na<sub>2</sub>-EDTA

Add ddH<sub>2</sub>O to 100 ml

Adjust pH to 8.0 with NaOH (use highly concentrated NaOH, Na<sub>2</sub>-EDTA will not dissolve, until the pH reaches 8.0)

Stored at room temperature

### **Acknowledgements**

BG11 medium is prepared according to Rippka *et al.* (1979). This work was supported by DFG grant to A. W. Wi2014/5-1.

### **References**

1. Dörrich, A. K., Mitschke, J., Siadat, O. and Wilde, A. (2014). [Deletion of the \*Synechocystis\* sp. PCC 6803 kaiAB1C1 gene cluster causes impaired cell growth under light-dark conditions.](#) *Microbiology* 160(Pt 11): 2538-2550.
2. Rippka, R., Deruelles, J., Waterbury, J. B., Herdman, M. and Stanier, R. Y. (1979). [Generic assignments, strain histories and properties of pure cultures of \*Cyanobacteria\*.](#) *J Gen Microbiol* (111): 1-61.

- 
3. Thiel, T., Bramble, J. and Rogers, S. (1989). [Optimum conditions for growth of \*Cyanobacteria\* on solid media.](#) *FEMS Microbiol Lett* 52(1-2): 27-31.