

Measurement of PI4P Levels in Intact Chloroplasts Isolated from *Arabidopsis thaliana*

Kumiko Okazaki^{1, 3*}, Shin-ya Miyagishima² and Hajime Wada¹

¹Department of Life Sciences, Graduate School of Arts and Sciences, The University of Tokyo, Tokyo, Japan; ²Department of Cell Genetics, National Institute of Genetics, Shizuoka, Japan;

³Present address: Graduate School of Bioscience and Biotechnology, Tokyo Institute of Technology, Yokohama, Japan

*For correspondence: okazaki.k.ac@m.titech.ac.jp

[Abstract] Phosphatidylinositol 4-phosphate (PI4P), a major species of phosphoinositides, modulates many fundamental cellular processes. We have recently revealed that PI4P plays an important role in chloroplast division as a negative regulator. Despite its importance in chloroplasts, the content of PI4P in chloroplasts is very low and it is difficult to measure PI4P levels. In this protocol, we describe a simple method that we have developed for measurement of low level of PI4P in chloroplasts. Intact chloroplasts were isolated by a basic method using Percoll gradient centrifugation and acidic lipids were extracted from the isolated chloroplasts. The extracted acidic lipids including PI4P were spotted onto the membrane strip, which had been pre-spotted with PI4P standards and other phosphoinositides as negative controls. PI4P in the spot of acidic lipids on the membrane was detected using a PI4P binding protein.

Materials and Reagents

A. Isolation of intact chloroplasts

1. Miracloth (Merck Millipore Corporation, Calbiochem[®], catalog number: 475855)
2. Gauze (Hakujiji Co., model: FC-gauze)
3. Seedlings of *Arabidopsis thaliana*
 Wild-type (Columbia-0) plants, *phosphatidylinositol 4-kinase (PI4K) α1* knockdown plants, *pi4kβ2-1* mutants and *PI4Kα1* knockdown plants of *pi4kβ2-1* mutants were grown for 4 days on MS agar plates and then transferred onto agar plates with or without 5 μM dexamethazone (DEX) and grown for 1 week. The down-regulation of *PI4Kα1* expression was induced by DEX treatments in *PI4Kα1* knockdown plants and *PI4Kα1* knockdown plants of *pi4kβ2-1* mutants. For treatments with inhibitors, wild-type 4-d-old seedlings were transferred onto agar plates with PI4K inhibitors, 200 μM wortmannin (WM) or 25 μM phenylarsine oxide (PAO), or a phosphatidylinositol 3-kinase (PI3K) inhibitor, 50 μM LY294002 (LY), or without inhibitors (DMSO) and grown for 3 days.
4. Percoll (GE Healthcare, Dharmacon, catalog number: 7-0891-01)
5. Bradford assay kit (Bio-Rad Laboratories, catalog number: 500-0006JA)
6. 1x protease inhibitor cocktail (Nakarai, catalog number: 03969-21)

Note: Currently, it is "Nacalai tesque, catalog number: 03969-21".

7. Sorbitol
8. HEPES-KOH (pH 7.5)
9. EDTA
10. Grinding buffer (see Recipes)
11. 80% Percoll or 40% Percoll (see Recipes)

B. Extraction of PI4P

1. 0.75 M Trichloroacetic acid (TCA) (Wako Pure Chemical Industries, Siyaku, catalog number: 203-04952)
2. 5% (w/v) TCA with 1 mM EDTA
3. Methanol (Wako Pure Chemical Industries, Siyaku, catalog number: 137-01823): Chloroform (Wako Pure Chemical Industries, Siyaku, catalog number: 038-02601) (2:1, v/v)
4. Methanol: chloroform: 12 N HCl (Wako Pure Chemical Industries, Siyaku, catalog number: 080-01066) (80:40:1, v/v/v)

C. Measurement of PI4P levels

1. PI(4)P Mass Strip Kit (Echelon Biosciences, catalog number: K-4000E)
2. Albumin from bovine serum (BSA), fatty acid free (Wako Pure Chemical Industries, Siyaku, catalog number: 017-15146)
3. HRP substrate solution (Thermo Fisher Scientific, Pierce, catalog number: NCI32132)
4. Phosphate (pH 7.4)
5. NaCl
6. 0.1% (w/v) Tween-20
7. Phosphate buffered saline (PBS) (see Recipes)
8. PBST (see Recipes)

Equipment

1. Homogenizer (Microtec Co., model: NS-51)
2. Centrifuge (TOMY SEIKO CO., model: MX-200 and GX-250)
3. Swing rotor (TS-7C)
4. Paintbrush
5. Pasteur pipette (Sansyo, Iwaki, catalog number: IK-PAS-9P)
6. Vacuum dryer (Centrifuge evaporator) (Shimadzu Corporation, model: SPE-200)
7. Sonicator (SHARP CORPORATION, model: UT-106)

Procedure

A. Isolation of intact chloroplasts

1. Sampling 0.5-1 g of seedlings of *Arabidopsis thaliana*. Whole seedlings for one-week-old plants or shoots for eleven-day-old plants were used.
2. Add 6 ml of grinding buffer and homogenize with homogenizer 3-5 times for 3 sec on ice.
3. Filtrate the homogenate through 2 layers of gauze and a layer of Miracloth (put gauzes on the Miracloth).
4. Centrifuge at $2,500 \times g$ for 15 min at 4 °C. Discard the supernatant and gently resuspend the pellet in 1.5 ml of grinding buffer with a paintbrush. Brush the pellet gently with paintbrush little by little to avoid breaking of intact chloroplasts. Note that pipetting the pellet breaks chloroplasts.
5. Add 1.5 ml of grinding buffer and 1 ml of 80% Percoll and gently mix. (Final concentration of Percoll is 20%.)
6. Prepare Percoll step gradients. Layering first 3 ml of 80% Percoll and then 5 ml of 40% Percoll in 15 ml centrifuge tubes. Put the tip of pipette on the surface and slowly lay 40% Percoll to avoid disturbance of the gradient. See Figure 1.
7. Lay 4 ml of crude chloroplast suspension (the results of step A5) on the top of Percoll gradients.
8. Centrifuge at $3,200 \times g$ for 30 min at 4 °C using a swing rotor. Set speed of acceleration and deceleration at minimum.
9. Collect a green band with a Pasteur pipette corresponding to intact chloroplasts at the interface between 40% and 80% to 2 ml tubes. Remove carefully the upper band (broken chloroplasts) and transfer lower band to 2 ml tubes with a Pasteur pipette. See Figure 1.
10. Add grinding buffer up to 2 ml and centrifuge at $700 \times g$ for 5 min at 4 °C.
11. Discard the supernatant and wash the pellet two times by resuspension in 1 ml of grinding solution and centrifugation at $700 \times g$ for 5 min at 4 °C.
12. Suspend the pellet (intact chloroplasts) in 0.3 ml of grinding buffer.
13. Take a part of the samples for measurement of protein concentration (step B9).

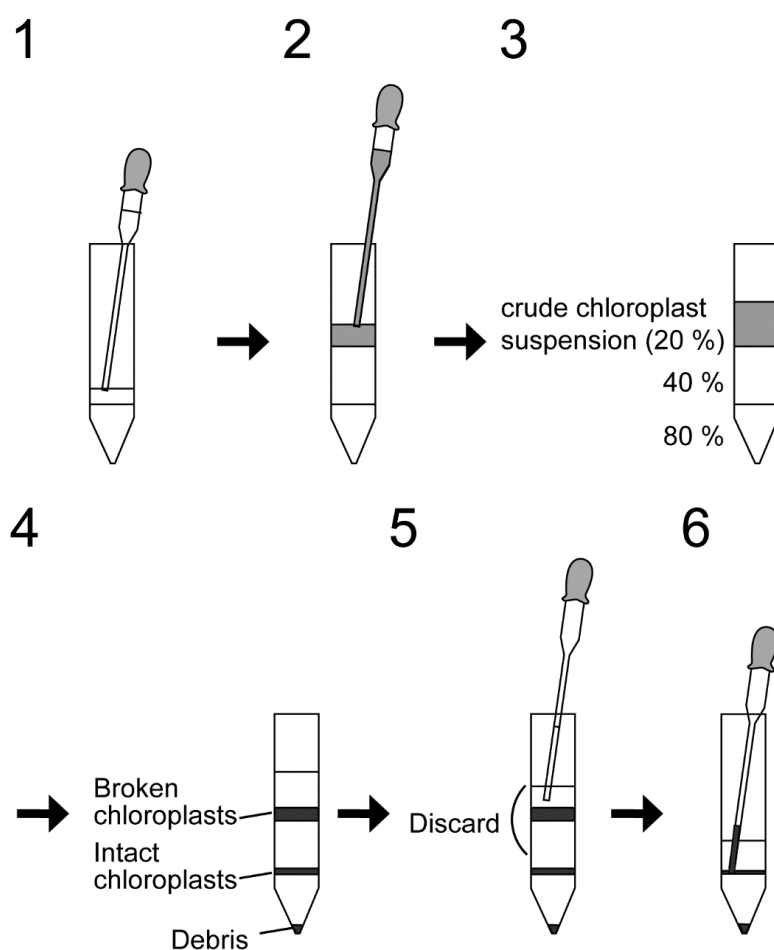


Figure 1. Isolation of intact chloroplasts with Percoll gradient centrifugation. 1. Lay 40% Percoll on 80% Percoll; 2. Lay crude chloroplast suspension on 40% Percoll; 3. The tube before centrifugation; 4. The tube after centrifugation; 5. Remove the upper band (broken chloroplasts); 6. Collect the lower band (intact chloroplasts).

B. Extraction of PI4P

1. Add 750 μ l of ice cold 0.75 M TCA to 250 μ l of intact chloroplast samples. Vortex and incubate on ice for 5 min.
2. Centrifuge at 700 x g for 5 min at 4 °C and discard the supernatant. Add 1 ml of 5% TCA with 1 mM EDTA to the pellet and vortex for 30 sec.
3. Centrifuge at 700 x g for 5 min at room temperature and discard the supernatant. Add 1 ml of methanol:chloroform (2:1) to the pellet for extraction of neutral lipids and vortex 3 times for 30 sec at an interval of 5 min at room temperature.
4. Centrifuge at 700 x g for 5 min at room temperature and discard the supernatant. Add 1 ml of methanol:chloroform (2:1) to the pellet for extraction of neutral lipids again and vortex for 1 min at room temperature.
5. Centrifuge at 700 x g for 5 min at room temperature and discard the supernatant. Add

750 μ l of methanol:chloroform:12 N HCl (80:40:1) to the pellet for extraction of acidic lipids and vortex 4 times for 30 sec at an interval of 5 min.

6. Add 250 μ l of chloroform and 450 μ l of 0.1 M HCl and vortex.
7. Centrifuge at 700 \times g for 5 min, collect the lower phase and transfer it to new tubes.
8. Dry samples in a vacuum dryer.
9. Measure protein concentration. Add 1-10 μ l of samples (step A13) into 1 ml of Bradford assay reagents (1:5 dilution). Incubate for 5 min at room temperature and measure absorbance at 595 nm. Measure absorbance of standard protein samples too. Make a standard curve using standard protein samples and calculate the PI4P concentration according to the reads. The concentration of protein is usually 0.3-0.8 mg/ml.

C. Measurement of PI4P levels

1. Add 15 μ l of chloroform: methanol: water (1:2:0.8, v/v/v) to the sample. Vortex for 30 sec and sonicate (100% output) in icy water bath for 5 min. Vortex again for 30 sec.
2. Centrifuge 15,000 \times g for 5 min and transfer 10 μ l of supernatant to new tube. Do not suck the pellet at the bottom of tube.
3. Spot samples onto a PI(4)P Strip. Spot 1 μ l of sample at a time. Before spotting sample again in the same area, dry the spot. After all samples are spotted, dry the strip for 30 min at room temperature. See Figure 2.
4. Incubate the membrane with blocking solution (3% fatty acid-free BSA in PBS) for 1 h at room temperature with gentle agitation.
5. Incubate the membrane with 5 ml of blocking solution containing 1 vial of PI(4)P Detector (2.5 μ g protein) for 1 h at room temperature with gentle agitation.
6. Discard the solution and wash the membrane with 10 ml of PBST three times at room temperature with gentle agitation for 5 min each.
7. Incubate the washed membrane with 5 ml of blocking solution containing 75 μ l of secondary detector for 1 h at room temperature with gentle agitation.
8. Discard the solution and wash the membrane three times with 10 ml of PBST at room temperature with gentle agitation for 5 min each.
9. Incubate the washed membrane with enough ECL solution to cover the membrane at room temperature for 5 min.
10. Wash the membrane with water three times to stop the reaction and dry the membrane.
11. Scan the membrane to record the obtained data. See Figure 2.

Representative data

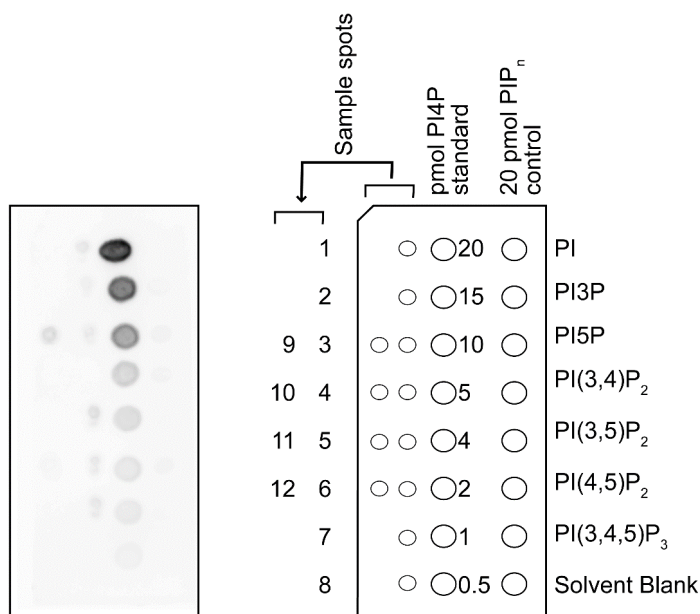


Figure 2. Representative results of PI4P detection assay. Samples extracted from intact chloroplasts were spotted left side of the PI(4)P Strip. PI4P standards (20, 15, 10, 5, 4, 2, 1, and 0.5 pmol of PI4P) and phosphoinositide controls were pre-spotted on the membrane. Lipids extracted from isolated chloroplasts of wild-type plants, *phosphatidylinositol 4-kinase (PI4K) α1* knockdown plants, *pi4kβ2-1* mutants and *PI4Kα1* knockdown plants of *pi4kβ2-1* mutants treated with/without 5 μM dexamethasone (DEX) were spotted onto the membrane. The down-regulation of *PI4Kα1* expression was induced by DEX treatments in *PI4Kα1* knockdown plants and *PI4Kα1* knockdown plants of *pi4kβ2-1* mutants. Lipids extracted from isolated chloroplasts of wild-type plants treated with/without PI4K inhibitors, 200 μM wortmannin (WM) or 25 μM phenylarsine oxide (PAO), or a phosphatidylinositol 3-kinase (PI3K) inhibitor, 50 μM LY294002 (LY), were also spotted. Sample spots are: 1. Wild-type plants without DEX; 2. Wild-type plants treated with DEX; 3. *PI4Kα1* knockdown plants without DEX; 4. *PI4Kα1* knockdown plants treated with DEX; 5. *pi4kβ2-1* mutants without DEX; 6. *pi4kβ2-1* mutants treated with DEX; 7. *PI4Kα1* knockdown plants of *pi4kβ2-1* mutants treated without DEX; 8. *PI4Kα1* knockdown plants of *pi4kβ2-1* mutants treated with DEX; 9. Wild-type plants without inhibitor; 10. Wild-type plants with WM; 11. Wild-type plants with PAO; 12. Wild-type plants with LY. The abbreviations used are: PI, phosphatidylinositol; PI3P, phosphatidylinositol 3-phosphate; PI5P, phosphatidylinositol 5-phosphate; PI(3, 4)P₂, phosphatidylinositol 3, 4-bisphosphate; PI(3, 5)P₂, phosphatidylinositol 3, 5-bisphosphate; PI(4, 5)P₂, phosphatidylinositol 4, 5-bisphosphate; PI(3, 4, 5)P₃, phosphatidylinositol 3, 4, 5-triphosphate.

Notes

1. The samples and Percoll gradient tubes should be kept on ice during procedure A.
2. Do not homogenize samples too much. Three to five times for 3 sec are usually enough. Longer homogenization results in lower yield of intact chloroplasts.
3. After step B3, do not keep samples on ice.
4. Intactness of chloroplasts can be assessed by phase contrast microscope. Intact chloroplasts are surrounded by bright halos of light. Broken chloroplasts look dark. See Figure 3 of Klinkenberg (2014).
5. Extraction of neutral lipids should be done twice for removing chlorophyll.
6. Do not suck the pellet at step C2. It contains left over proteins from extraction steps and makes hydrophobic film on the membrane. The hydrophobic film prevents binding of PI4P detector to PI4P. See Figure 3.



Figure 3. An example of a failed PI4P detection assay. Lipids were spotted onto the membrane together with the debris left over in the extraction step (the pellet at step C2). It prevented binding of PI(4)P Detector to PI4P and thus spots were detected as rings instead of circles.

Recipes

1. Grinding buffer
 - 0.33 M sorbitol
 - 30 mM HEPES-KOH (pH 7.5)
 - 2 mM EDTA
 - 1x protease inhibitor cocktail
 Autoclave without protease inhibitor cocktail
 Protease inhibitor cocktail is added before use.

2. 80% or 40% Percoll
 - 30 mM HEPES-KOH (pH 7.5)
 - 0.33 M sorbitol
 - 2 mM EDTA
 - Percoll (80% or 40%, v/v)
3. PBS
 - 10 mM phosphate (pH 7.4)
 - 150 mM NaCl
 - Sterilized by autoclave
4. PBST
 - PBS
 - 0.1% (w/v) Tween-20

Acknowledgments

This protocol was adapted from the original work (Okazaki *et al.*, 2015) to provide the detailed procedures. This work was supported by JSPS KAKENHI Grant Number 26840089.

References

1. Klinkenberg, J. (2014). [Extraction of chloroplast proteins from transiently transformed *Nicotiana benthamiana* leaves](#). *Bio-protocol* 4(18): e1238.
2. Nakanishi, H., Suzuki, K., Kabeya, Y. and Miyagishima, S. Y. (2009). [Plant-specific protein MCD1 determines the site of chloroplast division in concert with bacteria-derived MinD](#). *Curr Biol* 19(2): 151-156.
3. Okazaki, K., Miyagishima, S. Y. and Wada, H. (2015). [Phosphatidylinositol 4-phosphate negatively regulates chloroplast division in *Arabidopsis*](#). *Plant Cell* 27(3): 663-674.