

Preparation of Chloroplast Lipid Membrane and Lipid-Protein Interaction Assay

Kumiko Okazaki^{1, 3*}, Shin-ya Miyagishima² and Hajime Wada¹

¹Department of Life Sciences, Graduate School of Arts and Sciences, The University of Tokyo, Tokyo, Japan; ²Department of Cell Genetics, National Institute of Genetics, Shizuoka, Japan;

³Present address: Graduate School of Bioscience and Biotechnology, Tokyo Institute of Technology, Yokohama, Japan

*For correspondence: okazaki.k.ac@m.titech.ac.jp

[Abstract] Lipid-Protein interaction assay is a method to search lipids, which are bound with proteins *in vitro*. Since membranes that are spotted with chloroplast lipids such as monogalactosyldiacylglycerol (MGDG), digalactosyldiacylglycerol (DGDG), and sulfoquinovosyldiacylglycerol (SQDG) are not commercially available, we extracted these lipids from cyanobacterial cells and spotted them onto membranes. The prepared membranes could be used for lipid-protein interaction assay.

Materials and Reagents

A. Preparation of lipids for chloroplast lipid membrane

1. Amersham Protran Premium NC0.2 (GE Healthcare, Dharmacon, catalog number: AP-10600081)
2. TLC silica gel 60 plate (Merck Millipore Corporation, catalog number: 105721)
3. Fused silica capillary column, 0.25 mm x 50 m (Shinwa Chemical Industries, catalog number: HR-SS-10)
4. Cells of the cyanobacterium, *Synechocystis* sp. PCC 6803 (for lipid extraction)
5. Methanol (Wako Pure Chemical Industries, Siyaku, catalog number: 137-01823)
6. Chloroform (Wako Pure Chemical Industries, Siyaku, catalog number: 038-02601)
7. 0.01% primuline (Sigma-Aldrich, catalog number: 206865), in 80% acetone
8. Pentadecanoic acid (Sigma-Aldrich, catalog number: 91446-5G)
9. Methanolic HCl (Sigma-Aldrich, catalog number: 33050-U), dilute with methanol 3 times
10. Hexane (Wako Pure Chemical Industries, Siyaku, catalog number: 085-00411)

B. Lipid-Protein interaction assay

1. PIP Strips, Membrane Lipid Strips (Echelon Biosciences, catalog number: P-6001 and P-6002)
2. Hybri-Bag (Cosmo Bio, catalog number: S-1001)
3. Albumin from bovine serum (BSA), fatty acid free (Wako Pure Chemical Industries, Siyaku, catalog number: 017-15146)

4. Purified tagged proteins
5. An anti-penta-His mouse monoclonal antibody (QIAGEN, catalog number: 34660) and an anti-GST mouse monoclonal antibody (Sigma-Aldrich, catalog number: G1160)
6. Horse radish peroxidase (HRP)-conjugated goat anti-mouse secondary antibody (Thermo Fisher Scientific, catalog number: 31432)
7. HRP substrate solution (Thermo Fisher Scientific, Pierce, catalog number: NCI32132)
8. Phosphate (pH 7.4)
9. NaCl
10. 0.1 % Tween-20
11. Phosphate buffered saline (PBS) (see Recipes)
12. PBST (see Recipes)

Equipment

1. Glass tubes (Iwaki, catalog number: 71-088-004 and 71-063-006)
2. Pasteur pipette (Sansyo, Iwaki, catalog number: IK-PAS-9P)
3. Spectrophotometer (Scintec Instruments, Pharmacia Biotech, model: Ultrospec 2000)
4. Vortex mixer
5. Centrifuge (TOMY SEIKO CO., model: LC-120)
6. Vacuum dryer (Shimadzu Corporation, model: SPE-200 centrifuge evaporator)
7. TLC developing chamber (AS ONE, catalog number: CK-0544-060)
8. Hair dryer (commercially available)
9. UV illuminator
10. Razor blade
11. Heat block
12. Gas chromatograph (Shimadzu Corporation, model: GC-18A)
13. Fused silica capillary column, 0.25 mm x 50 m (Shinwa Chemical Industries, catalog number: HR-SS-10)
14. Data processor (Shimadzu Corporation, model: C-R2AX)

Procedure

- A. Preparation of chloroplast lipid membrane
 1. Collect cyanobacterial cells by centrifugation and suspend the collected cells in 2 ml of 0.9% KCl. 3×10^9 cells of *Synechocystis* PCC 6803 are sufficient [Number of cells were determined by spectrophotometer. Optical density at 730 nm (OD_{730}) = 1.0 corresponds to 1×10^8 cells /ml for *Synechocystis* PCC 6803].
 2. Add 7.5 ml of methanol:chloroform (2:1, v/v) and vortex for 1 min.
 3. Add 2.5 ml of chloroform and 2.5 ml of H₂O, and vortex.

4. Centrifuge at 700 x g for 10 min. Collect lower phase, put it into clean tube and dry in a vacuum dryer. Measure the weight of the tube before use and after drying lipids. Calculate the yield of lipids.
5. Dissolve dried lipids in chloroform at 10 mg/ml.
6. Add developing solvent into a TLC developing chamber, close and leave it for 30 min as shown in Figure 1. Developing solvent is acetone:benzene:ultrapure water = 91:30:8 (v/v/v). Spot the chloroform solution containing 1-3 mg of lipids onto a TLC plate and place the plate into the TLC developing chamber.

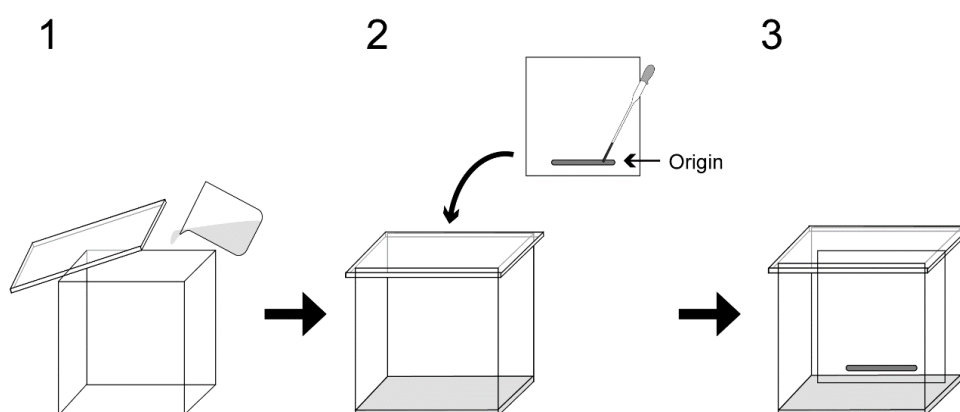


Figure 1. Separation of lipids on TLC plate. 1. Add developing solvent into a TLC developing chamber; 2. Spot lipid solution onto a TLC plate and place the plate into the chamber; 3. The TLC plate is developed in the chamber to separate lipids.

7. When the solvent front has reached 1 cm from the top of the plate, remove the plate from the chamber and dry with cool wind provided by hair dryer.
8. Spray 0.01% primuline in 80% acetone onto the plate. Lipids are viewed as spots under long wavelength of UV light as shown in Figure 2. Mark the corresponding position of each lipid and scrape the silica gel in the each spot from the TLC plate that is grass plate coated with silica gel. The silica gel can be scraped from the plate with a razor blade.

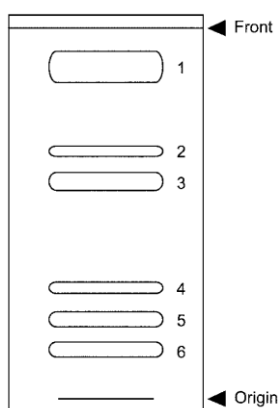


Figure 2. Thin-layer chromatogram of cyanobacterial lipids visualized by UV.

Total lipids extracted from cyanobacteria were developed in acetone:benzene:ultrapure water=91:30:8. 1. Pigments; 2. Monoglucosyldiacylglycerol; 3. Monogalactosyldiacylglycerol (MGDG); 4. Digalactosyldiacylglycerol (DGDG); 5. Sulfoquinovosyldiacylglycerol (SQDG); 6. Phosphatidylglycerol (PG).

9. Add 1 ml of chloroform to the silica gel scraped from the plate and vortex. Spin down the silica gel and transfer supernatant to a new tube. Repeat this extraction 3 times and evaporate chloroform in a vacuum dryer. Dissolve the obtained lipids in 100 μ l of chloroform. The following 10 to 12 steps are required to determine the yield of obtained lipids.
10. Take 5 μ l of the chloroform solution containing each lipid to a new tube with a screw cap and add 50 nmol of pentadecanoic acid as an internal standard and 2.5 ml of 3% (w/v) methanolic HCl. Close the tube with the screw cap and incubate at 80 $^{\circ}$ C for 3 h. During the incubation fatty acids bound to lipids are derivatized to fatty acid methyl esters.
11. After cooling the tube to room temperature, add 2.5 ml of hexane and vortex vigorously. Fatty acid methyl esters are extracted with hexane. Leave to stand for 5 min and collect upper phase (hexane phase), put it into a clean tube and evaporate hexane in a vacuum dryer.
12. Dissolve the obtained fatty acid methyl esters in 50 μ l of hexane and analyze them with a gas chromatograph. Temperatures for the injector, column and detector chambers were set at 250 $^{\circ}$ C, 170 $^{\circ}$ C and 250 $^{\circ}$ C, respectively. Determine lipid content by calculating total amount of fatty acid methyl esters based on the peak areas of fatty acid methyl esters in the obtained chromatograms. Pentadecanoic acid is not contained in cyanobacterial lipids, thus it can be used as an internal standard.
13. Dilute lipids (the obtained lipids at the step A9) with chloroform to 0.1, 1 and 10 mg/ml.
14. Spot the lipids corresponding to 0.1, 1 and 10 mg/ml in 1 μ l of chloroform onto

Amersham Protran Premium membrane and dry completely.

B. Lipid-Protein interaction assay

1. Incubate the membrane with 5 ml of blocking solution (3% fatty acid-free BSA in PBS) for 1 h at room temperature with gentle agitation.
2. Discard the blocking solution and incubate the membrane with 2 ml of blocking solution containing 0.5 µg/ml tagged-protein (His-tagged or GST-tagged protein) overnight at 4 °C. Put the membrane and solution into a Hybri-Bag and seal it.
3. Discard the protein solution and wash the membrane with 10 ml of PBST three times at room temperature with gentle agitation for 10 min each.
4. Discard the solution and incubate the membrane with 2 ml of blocking solution containing primary antibody (anti-penta-His mouse antibody 1:2,000 or anti-GST mouse monoclonal antibody 1:2,000) for 1 h at room temperature with gentle agitation. Put the membrane and solution into a Hybri-Bag and seal it.
5. Discard the solution and wash the membrane with 10 ml of PBST three times at room temperature with gentle agitation for 10 min each.
6. Discard the solution and incubate the membrane with 2 ml of blocking solution containing secondary antibody (HRP-conjugated goat anti-mouse antibody 1:20,000) for 1 h at room temperature with gentle agitation. Put the membrane and solution into a Hybri-Bag and seal it.
7. Discard the solution and wash the membrane with 10 ml of PBST three times at room temperature with gentle agitation for 10 min each.
8. Incubate the membrane with enough HRP substrate solution to cover the membrane (2 to 3 ml) at room temperature for 5 min.
9. Wash the membrane with ultrapure water to stop the reaction and dry the membrane.
10. Scan the membrane to record the obtained data.

Representative data

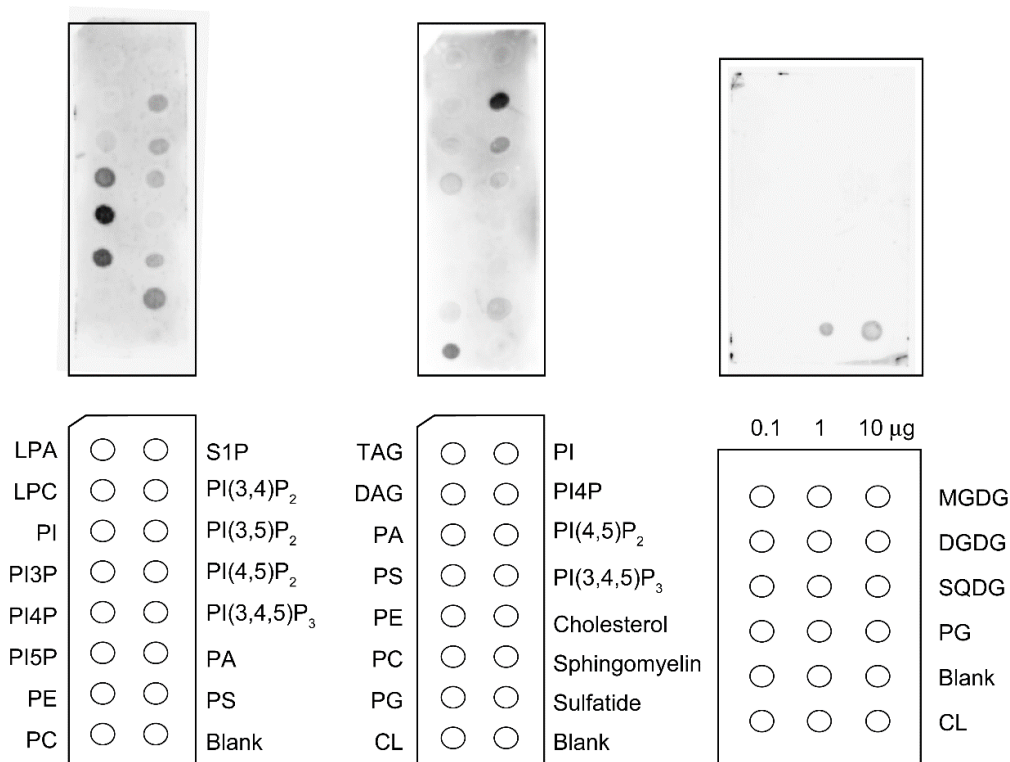


Figure 3. Representative results of Lipid-Protein Interaction Assay. Each spot on the left (PIP Strip) and middle membrane (Membrane Lipid Strip) contained 100 pmol of lipids. The right membrane contained chloroplast lipids, MGDG, DGDG, SQDG, and PG, and also contained cardiolipin (CL) as the control. The membranes were incubated with 0.5 µg/mL GST-PDV2. The abbreviations used are: LPA, lysophosphatidic acid; LPC, lysophosphatidylcholine; PI, phosphatidylinositol; PI3P, phosphatidylinositol 3-phosphate; PI4P, phosphatidylinositol 4-phosphate; PI5P, phosphatidylinositol 5-phosphate; PE, phosphatidylethanolamine; PC, phosphatidylcholine; PI(3, 4)P₂, phosphatidylinositol 3, 4-bisphosphate; PI(3, 5)P₂, phosphatidylinositol 3, 5-bisphosphate; PI(4, 5)P₂, phosphatidylinositol 4, 5-bisphosphate; S1P, sphingosine 1-phosphate; PI(3, 4, 5)P₃, phosphatidylinositol 3, 4, 5-triphosphate; PA, phosphatidic acid; PS, phosphatidylserine; TAG, triacylglycerol; DAG, 1, 2-diacylglycerol. The image was reprinted from Reference 2 with permission from *The Plant Cell*; copyright American Society of Plant Biologists (www.plantcell.org).

Recipes

1. PBS
 - 10 mM phosphate (pH 7.4)
 - 150 mM NaCl

2. PBST
- PBS
- 0.1 % Tween-20

Acknowledgments

This protocol was adapted from the original work (Okazaki *et al.*, 2015) to provide the detailed procedures. This work was supported by JSPS KAKENHI Grant Number 26840089.

References

1. Bligh, E. G. and Dyer, W. J. (1959). [A rapid method of total lipid extraction and purification](#). *Can J Biochem Physiol* 37(8): 911-917.
2. Okazaki, K., Miyagishima, S. Y. and Wada, H. (2015). [Phosphatidylinositol 4-phosphate negatively regulates chloroplast division in *Arabidopsis*](#). *Plant Cell* 27(3): 663-674.
3. Sato, N. and Murata, N. (1988). [Membrane lipids](#). *Meth Enzymol* 167: 251-259.